

Focus On ALIGNER SYSTEMS






3M Oral Care
(800) 423-4588
3m.com/clarityaligners

**Accusmile Digital
Lab Service Group**
(800) 721-4940
accusmile.de/en/
orthodontists/

**Align Technology
Inc**
(888) 822-5446
invisalign.com/partner

**Allesee Orthodontic
Appliances**
(800) 262-5221
aoaaccess.com

DynaFlex
(800) 489-4020
dynaflex.com

					
Product name	3M Clarity Aligners	Accusmile Aligners	Invisalign Clear Aligners	Simpli5	EZ-Align Clear Aligners
What materials is the product made of?	Proprietary	TRACK thermoforming foils from Forestadent	Invisalign SmartTrack material—a proprietary multi-layered polymer designed specifically for use with the Invisalign system.	Proprietary	Made from Glacier, an ultra premium proprietary thermoformed material that focuses on patient comfort and strength and durability.
Made in lab or in-house	In-house	Both	In-house	Lab	Both
Software interface	The 3M Oral Care Portal is a HIPAA compliant, open, web-based platform allowing practitioners access to create, refine, and update treatment plans at anytime from anywhere. The portal is compatible with any market-leading intraoral scanner, giving practitioners the control and flexibility needed when treatment planning for patients.	OnyxCeph ³	ClinCheck Pro with 3D Controls	N/A	Proprietary new Smile-Share software uses the newest technology for case management and tooth movement changes. Offers direct to manufacture after doctor approval of treatment plan or direct changes made by the doctor.
Orthodontist's cost per case	Contact company. Varies by annual volume. Discounts available.	Varies by product and volume	Varies by product	Contact company. Discounted pricing on single arch cases.	Contact company. Varies by system.
Special features	The 3M Clarity Aligners System delivers custom aesthetic treatment options that showcase an orthodontist's expertise. The 3M Digital Bonding system is designed to provide added efficiency and accuracy throughout the bonding process with the added convenience of being able to create customized combination treatment plans in the 3M Oral Care Portal.	The Accusmile service from Forestadent allows orthodontic practices the need-based digital planning and fabrication of aligners. From digitizing the models, treatment planning, 3D printing of the arches to deep-drawing the aligners — you are always in control of the process.	Weekly aligner changes may reduce treatment time up to 50% compared to 2-week changes. Precision cuts for use of elastics—designed to eliminate need to manually cut aligners. Depending on package, additional aligners may be available at no additional charge for up to 5 years maximum per case.	Designed to be simple and reliable, Simpli5 has been a dependable option for Express aligner treatment for anyone looking to perfect a beautiful smile. No approver required.	The product was recently updated in its software and material to offer the highest quality product. SmileShare software offers the ability for the doctor to fully customize each treatment plan and the overall experience. Glacier material utilizes undercut to fully engage the teeth for maximum results. The cut of the aligners also offers the most support during the tooth movement process.

Information in this guide based on data submitted by product manufacturers. Orthodontic Products strives for accuracy in all data but cannot be held responsible for claims made by manufacturers.

**Henry Schein
Orthodontics**
(800) 344-3347
revealclearaligners.
com/orthodontists/


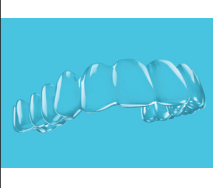




Ormco Corp
(800) 854-1741
ormco.com/spark

**SML - Space Main-
tainers Laboratories**
(800) 423-3270
smlglobal.com

**Tortorici Ortho/
Clear Fusion**
(631) 642-1900
clearfusionortho.com

TP Orthodontics Inc
(800) 348-8856
tportho.com

**Tru-Tain
Orthodontic &
Dental Supply**
(800) 843-0904
tru-tain.com

					
Reveal Plus Clear Aligner System	Spark Clear Aligners	Spring Retainer	Clear Fusion Hybrid Retainer	Refine Complete Aligner System	Tru-Tain DX
Proprietary ClearWear material for ultimate clarity and fit, reducing the need for attachments.	TruGEN and TruGEN XR offer greater sustained force retention. Designed for more effective tooth movement than the leading aligner brand.	The Spring Hawley Retainer is useful in correcting minor rotations and crowding of up to 1 mm to 3 mm, cuspid to cuspid, on the upper or lower arch.	Clear thermal reinforced with orthodontic acrylic	Thermoplastic polymer	Hybrid blend copolyester plastics
Lab	In-house	Lab	Lab	Lab	Both
The intuitive Digital Dental Exchange (DDX) Reveal Portal streamlines the case submission and treatment planning process.	Spark's 3D Approver software is designed to give doctors more control and flexibility for more efficient and predictable treatment planning.	SML can provide IPR analysis; extra fees apply.	N/A	TPRX Portal (tportho.com/custom-lab-services/)	N/A
Contact company	Volume-based discounts with the base at \$1,649.	\$101.20	\$65/arch	Contact company. Treatment times vary.	N/A
Reveal Plus Clear Aligners are designed uniquely clear, provide a precision fit for fewer needed attachments, and offer numerous workflow and treatment efficiencies. Reveal Plus incorporates the paradigm changing SAGITTAL FIRST Philosophy with Carriere MOTION 3D CLEAR Appliances to achieve ultimate aesthetics and efficiency.	Spark Clear Aligners introduced new features not previously available with clear aligner therapy. These include TruGEN XR, a new material option for finishing and refinements, the ability to add Posterior Bite Turbos during treatment planning, and CBCT integration, enabling orthodontists to import and visualize CBCT patient data within Spark's software.	Space necessary to correct this crowding is gained with IPR in the anterior region only. The appliance is fabricated to this corrected anterior position. On the day of delivery, complete the necessary IPR, then place the appliance. The spring action of the cuspid wires will direct a light labial-lingual force to align the teeth. Fabricated to the "ideal" setting, the appliance may not seat well in the patient's mouth—until the teeth adapt to their new position.	Rigid clear retainer	The Refine Complete Aligner System addresses end-of-treatment challenges, such as posterior open bite commonly seen in aligner patients. Treatment options are designed for flexibility, and each case includes your choice of a custom or semi-custom finishing appliance. No advanced training is required—cases can be easily submitted through TPRX.	Durable retainer and aligner plastic